



**Functional Protection**  
Without Compromise



# DONJOY® Functional Protection Without Compromise

One device for multiple stages, from recovery to rehabilitation through return to activity.

## The DonJoy Functional Wrist Brace

is the first functional post-injury wrist brace offering the unprecedented combination of rigid support and range-of-motion control providing performance and protection for any activity—without compromise.



### Ability to Control Range-of-Motion

#### Variable Extension Control

The patented articulating hinge system allows for controlled and anatomically accurate movement of the wrist. The removable extension stops control motion at **0°, 20°, 40°** and **60°**, based on individual needs progressing from immobilization through rehabilitation and back to work or sport.

### Light, Comfortable Fit

#### Cool Material

High-grade, open-cell liner provides a cool, comfortable fit. Material is washable and waterproof.

### Low-Profile Design

Designed to fit underneath sporting equipment, the brace is comfortable throughout the day.



### Radial/Ulnar Control

An optional radial ulnar deviation stop can be inserted to control side-to-side motion as needed.

### Open Palm for Easy Gripping


A “top-loaded” dorsal frame leaves the palm open, allowing performance of activities that require a hand-gripping motion.

### Removable Stay

A moldable palmer stay is inserted when additional stability is needed and is removed to allow mobility.



**CXT:** Carbon fiber and Kevlar®, for maximum strength

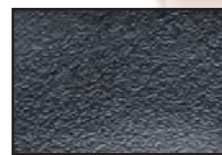

**SXT:** Injection molded from composite material

**Custom:** Options are available in CXT model only, and include a variety of colors



**SPEED** A speed wrap for quick application—lowest profile design

**SPORT** A forearm strap pulls through the brace for support and stability

**LACER** Easy-pull lacing system that connects to palmer stay—best for longer periods of wear

# 3 CONFIGURATIONS



**DONJOY**

1430 Decision Street Vista, CA 92081-8553 U.S.A.  
T 800.336.6569 D 760.727.1280 F 800.936.6569  
[www.donjoy.com](http://www.donjoy.com)



[DJOglobal.com](http://DJOglobal.com)

DJO proudly sponsors:



Kevlar is a registered trademark of E.I. du Pont de Nemours and Company.

© 2010 DJO, LLC 00-0889 Rev A